


Forever Medical Autoclave

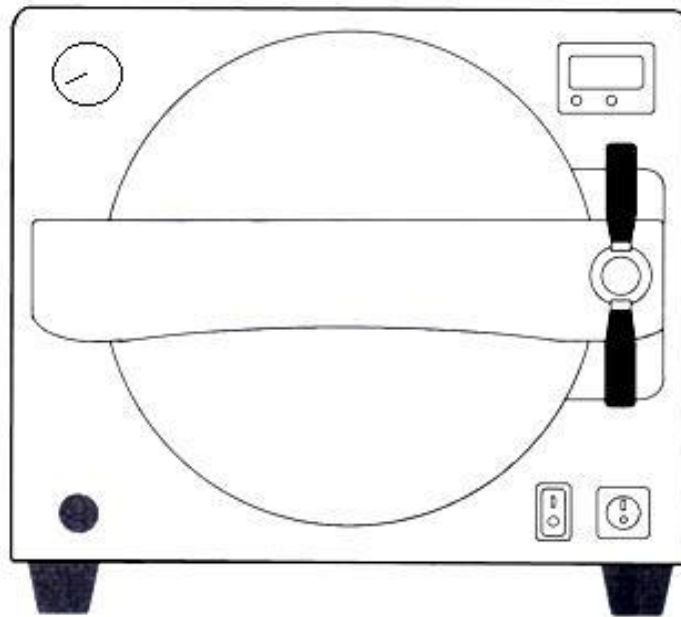

Model No:YJ-AN18

USER`S MANUAL

Thank you very much for purchasing our product.

Please read and understand all the instructions before using this machine and keep this manual for your timely reference.

SPECIFICATION

- CHAMBER SIZE : stainless steel #304 ; 18Liter ; (Ø)25×(D)35 CM

- MACHINE SIZE : (L)58×(W)39×(H)36 C

- POWER : AC220V/900W/4A

- TEMP/PRESSURE : 121℃/0.12Mpa ; 134℃/0.22Mpa

- STERILIZATION TIME : 25Min / 6Min

- SAFETY : Heating-over protection ; safety valve

- TEMP CONTROL DESIGN : chip control

- OPTIONAL ACCESSORIES : tray


basket × 1pc ; tray holder × 1pc

Sterilization tray × 3pcs ; steam exhaust pipe × 1pc


P.1

NAME OF EACH PART

- | | |
|--------------------------------|--------------------------|
| 1. Temp display window | 2. 134°C light |
| 3. Error indicator 134°C light | 4. Door handle |
| 5. Power switch | 6. Temp selector switch |
| 7. Pressure gauge | 8. Steam discharge valve |
| 9. Steam vent | 10. Safety valve |

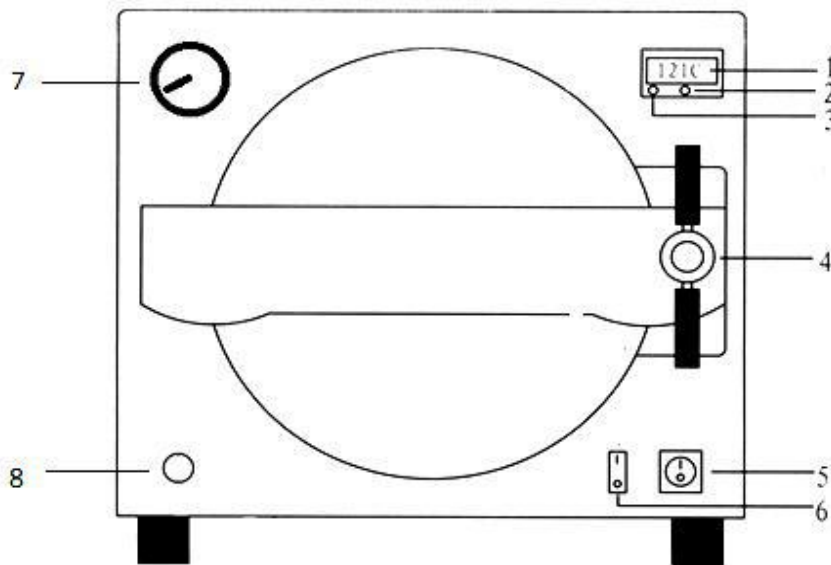


Back of mahchine:

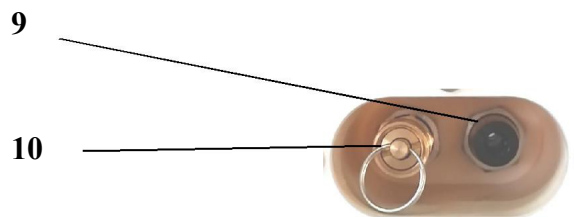

INSTALLATION

The machine is put in the absence of inclined plane on the desktop.

The machine please the postiton to be away from wall 5cm.

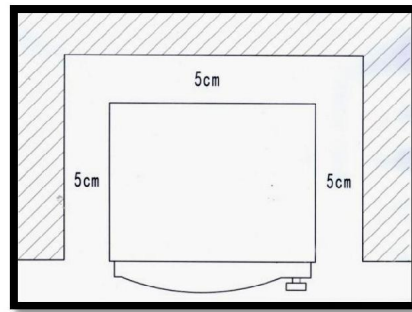

P.2

USE STEP

Step 1: into about 800 ml of distilled water disinfection chamber.

Step 2: put the disinfection equipment in the tray.


Step 3: turn right for close the steam discharge valve.

Step 4: choose your need sterilization temperature. (121 °C or 134 °C)

Step 5: open the power switch.


Step 6: close the door, the door handle to turn right.

(Rotate to a safe location, start the power switch. Get to work).

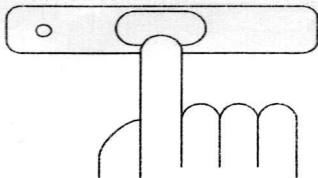

CAUTION

1. After the completion of disinfection。 Please to wait for 10 minutes。
Before open the DOOR(Open the discharge valve of steam turn left) ◦

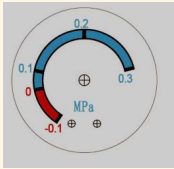

2. pressure Gauge in **red area**. Do not rotate handle.

The first to open the steam discharge valve .

Gauge pressure to **zero position**, such as in open the door

P.4

TROUBLESHOOING

When disinfection work has not completed,the buzzer makes a sound and the Error indicator twinkle.

Please to check: THE FIRST ,turns left the Steam drain valve. emissions steam and pressure.Confirm the fault codes temperature display window

E1 & E4: CHAMBER dry (lack of water)

off the power switch, waiting for 10 minutes. Let CHAMBER temperature drop.Open the door and add enough water.

E2: Temperature sensor fault

E3: The supply of voltage instability (restart the machine again)